


Bowie Dick Test Sheet

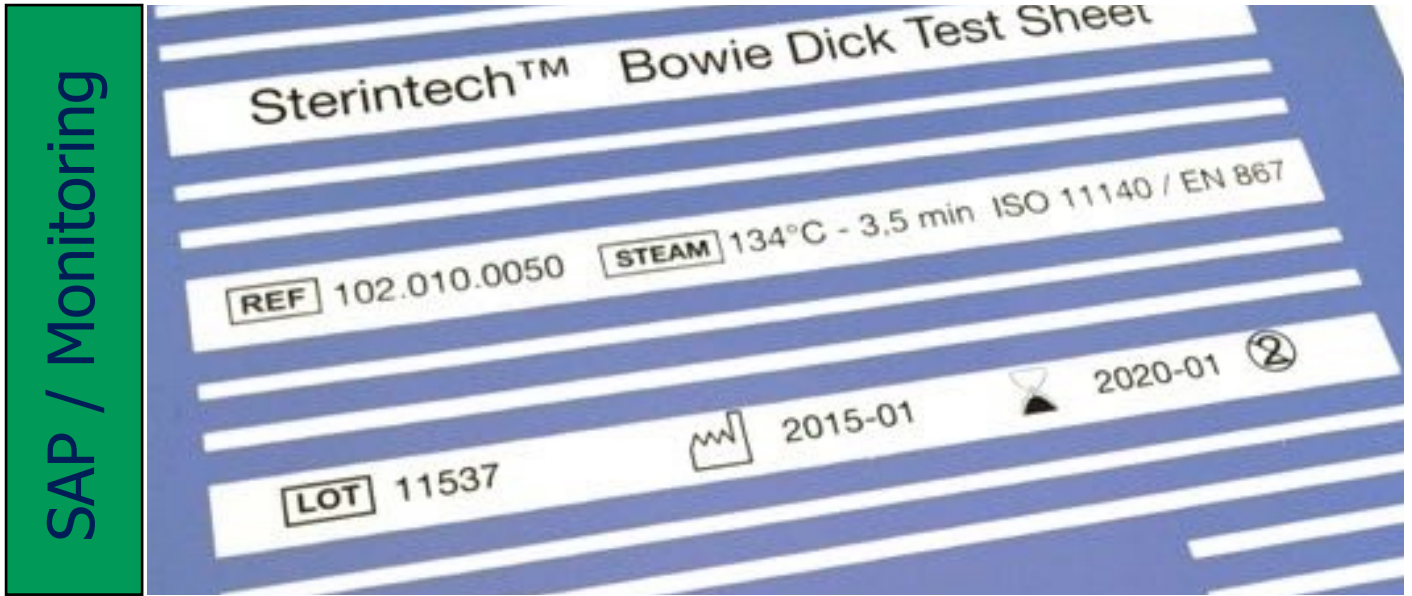

Once a day every day

REF: 102.010.0050 Bowie Dick Test Sheet A4 Unit: 50 pieces
REF: 102.011.0050 Bowie Dick Test Sheet A5 Unit: 50 pieces

Used for daily monitoring of pre-vacuum steam sterilizer. Based upon the Bowie Dick test result the steam sterilizer is released for usage during the day.

Characteristics:

- ◆ For 134 °C – 3,5 min / 121 °C – 15 min
- ◆ To be packed into pack as per standard / instructions
- ◆ Fulfill requirements of EN 285
- ◆ Conform to ISO 11140 - 1:2014 Type 2
- ◆ Dimension: DIN A4 size
- ◆ Clear Blue to Pink color change
- ◆ Easy to interpret / Easy to document
- ◆ Non-Toxic / Lead Free

The Bowie Dick test detect failures like air leaks, insufficient vacuum, inadequate steam penetration.